
nuance audio

Nuance Audio Remote Plus

Remote Plus guided onboarding
Setup Flow

The Remote Plus is used to activate and set up your eyewear. It guides you through onboarding, lets you choose your preferred language, and helps you through a quick personal calibration so everything feels just right for you.

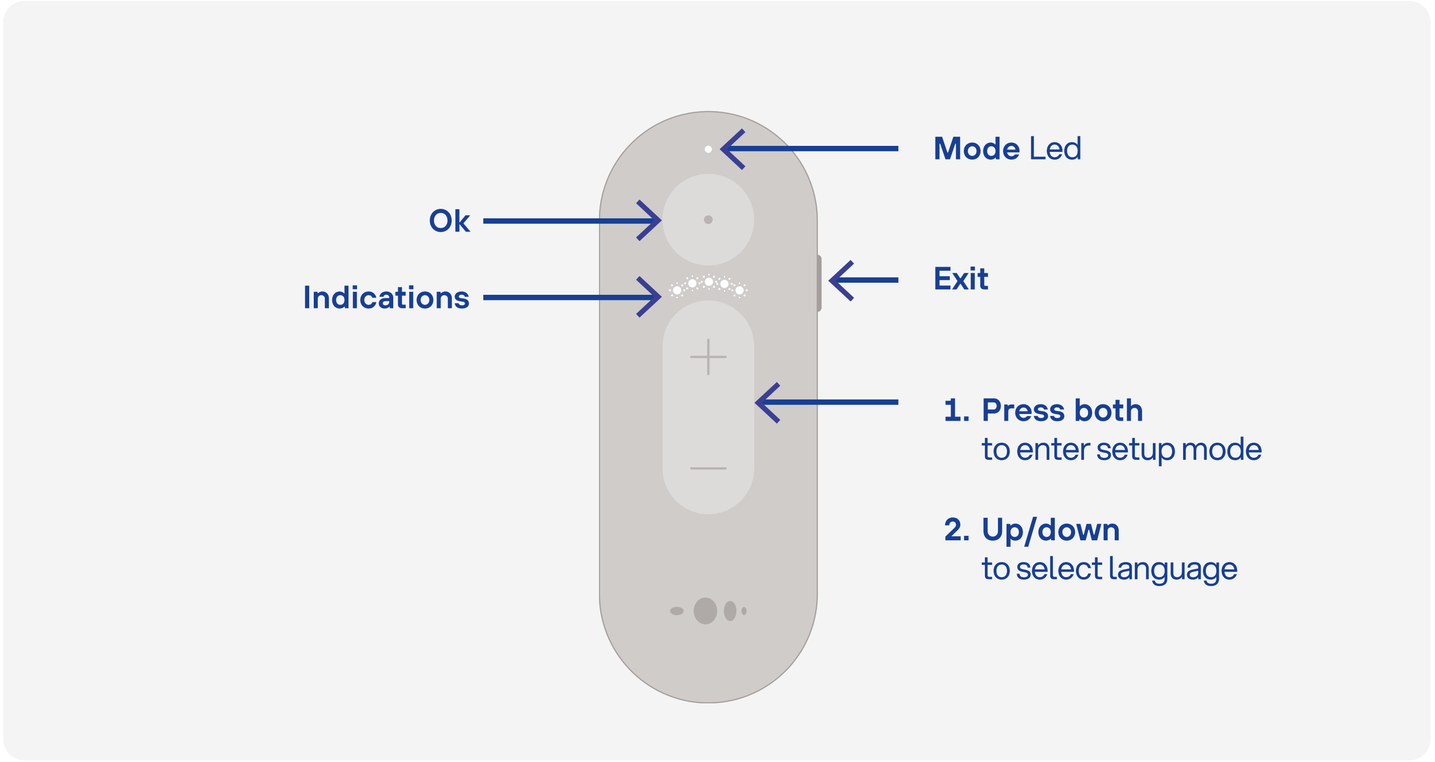

Remote Plus guided onboarding
Setup Flow

Get ready

1. Glasses and Remote Plus are paired
2. Glasses are on and worn
3. Quiet room

Pairing options
1. *Pre-paired hard bundle*
2. *Store Associate pairs in store*
3. *User pairs at home*

Enter setup mode

1. Press and hold **+/-** buttons

Glasses measure environmental noise to proceed

Select language

1. Hear languages
2. Scroll with **+/-**
3. Confirm with top button

Audio prompts are used for calibration feedback

Calibrate


1. Glasses:
   “Repeat after me: 1, 2, 3, 4”
2. User: “1, 2, 3, 4”
3. Continue...

- *Audio prompts are used for calibration feedback*
- *User voice is used for own-voice calibration*


Done

1. Setup complete once enough data has been collected

Remote Plus user interface
Frame Pairing and Set Up

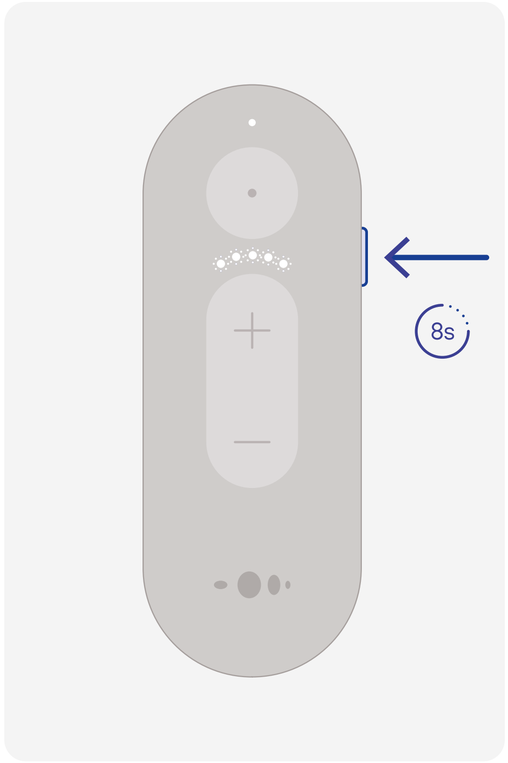


Pairing
Remove the battery tab from the Remote before proceeding. When the Glasses are on, press and hold the "**+**" or "**-**" volume button on the frame until blue LED blinks. Then, **press side button** and hold for 8 seconds.

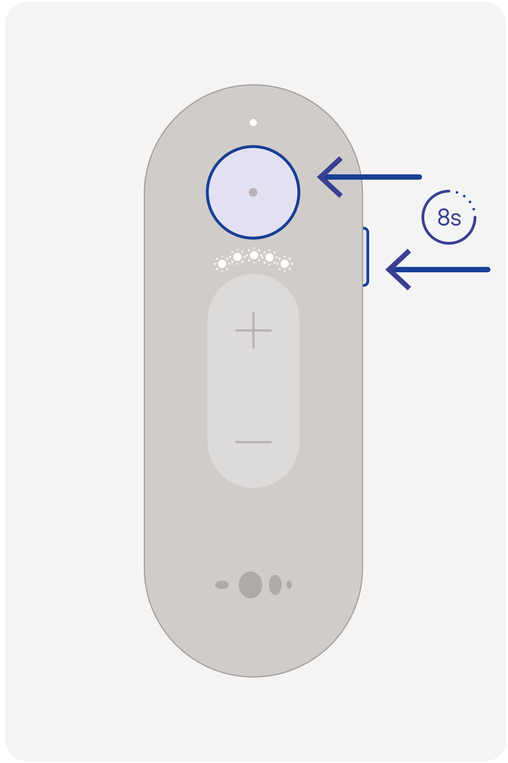


Factory Reset
Press top button and side button. Hold both for 8 seconds.


Remote Plus user interface
Frame Pairing and Set Up

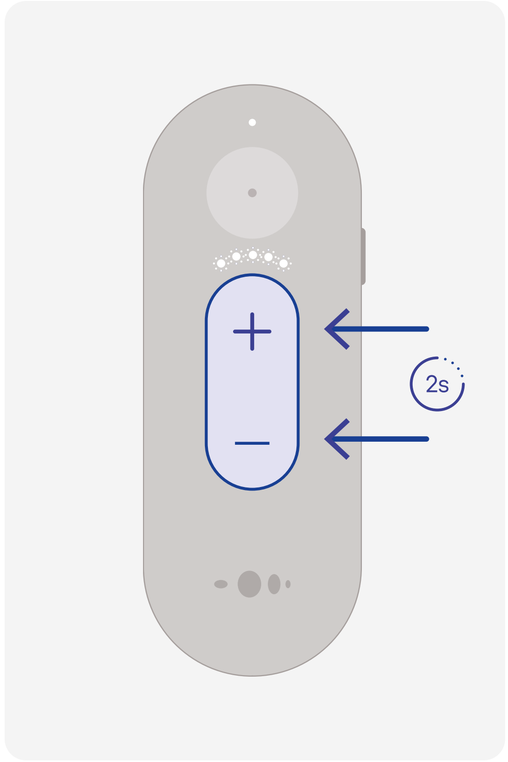

Enter setup mode
Press and hold both (+) and (-) for 2 seconds.

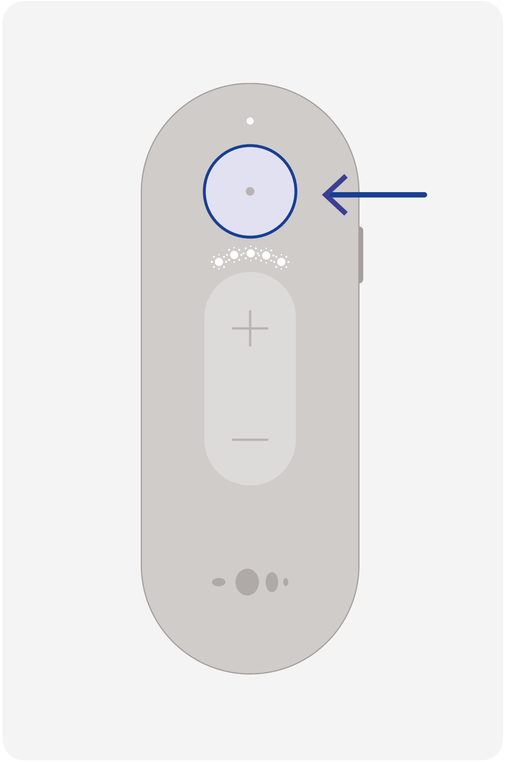


OK
Press top button

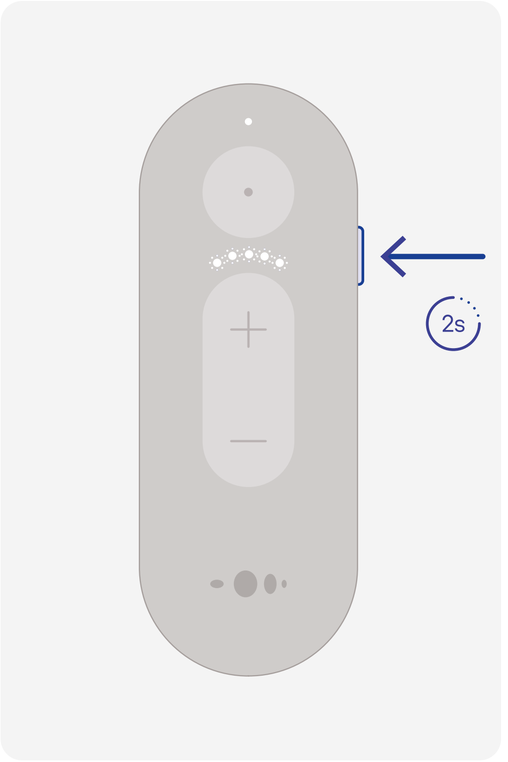


Exit setup
Press side button and hold for 2 seconds.

Remote Plus user interface
Frame Setup

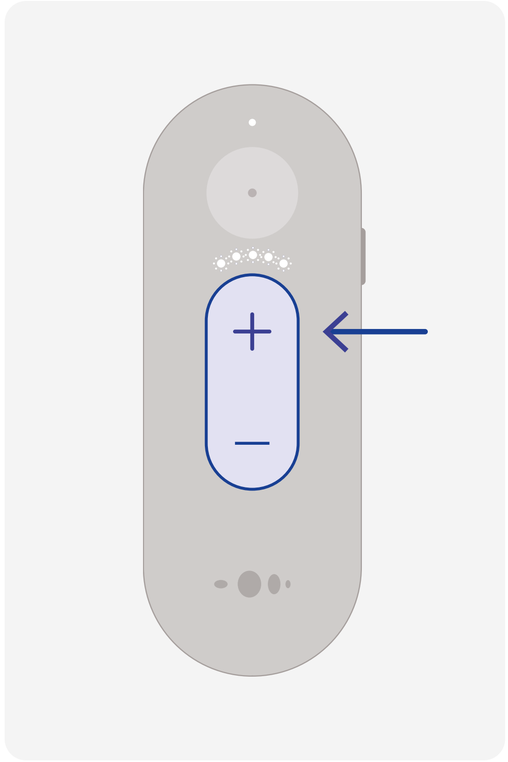


Next language
Press (**+**)

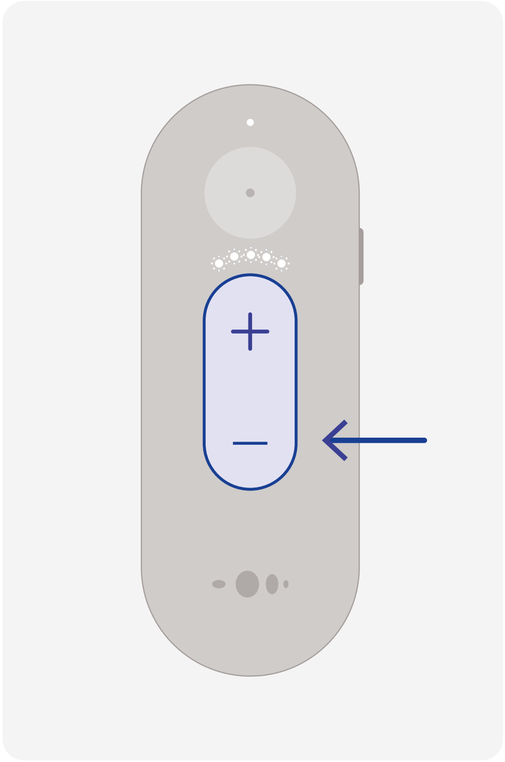


Previous language
Press (**-**)

Remote Plus interface
Volume

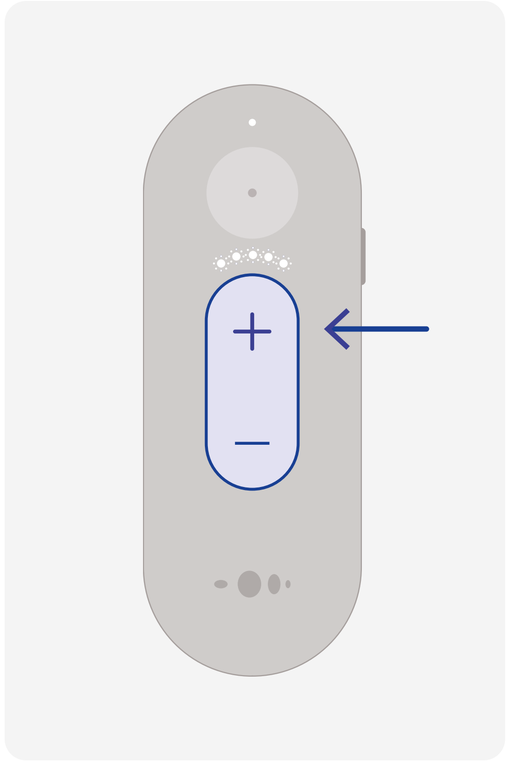


Volume up
Press (+)

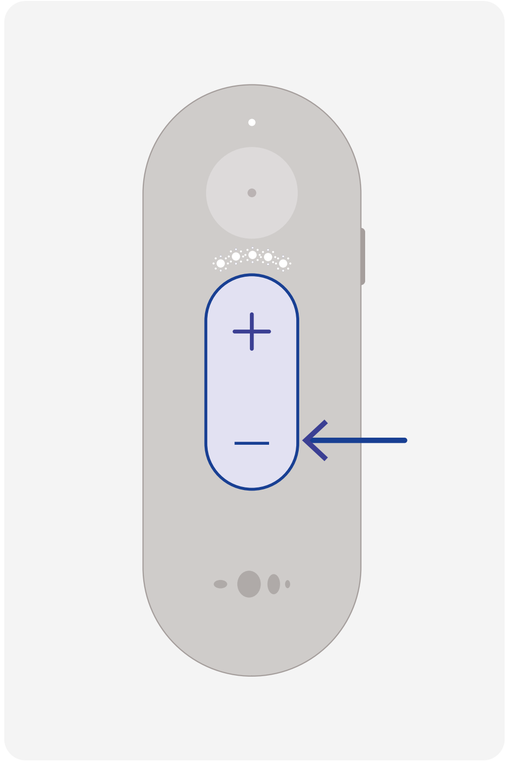


Volume down
Press (-)

Remote Plus user interface
Volume

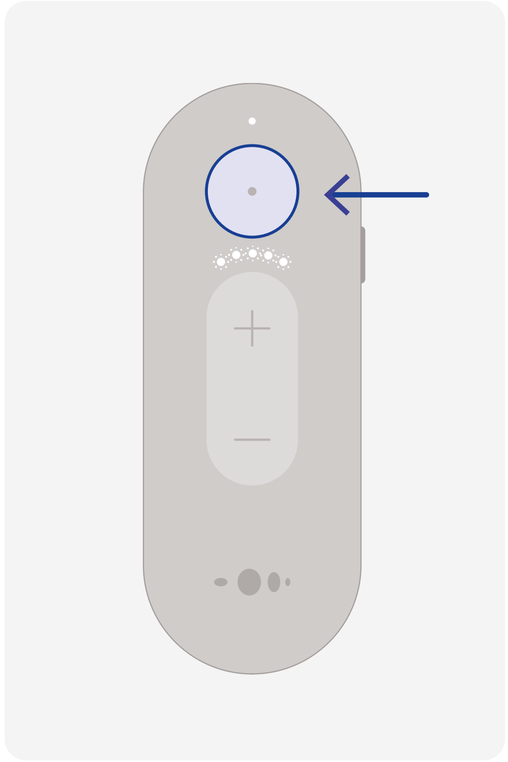


Switch audio mode
Press top button

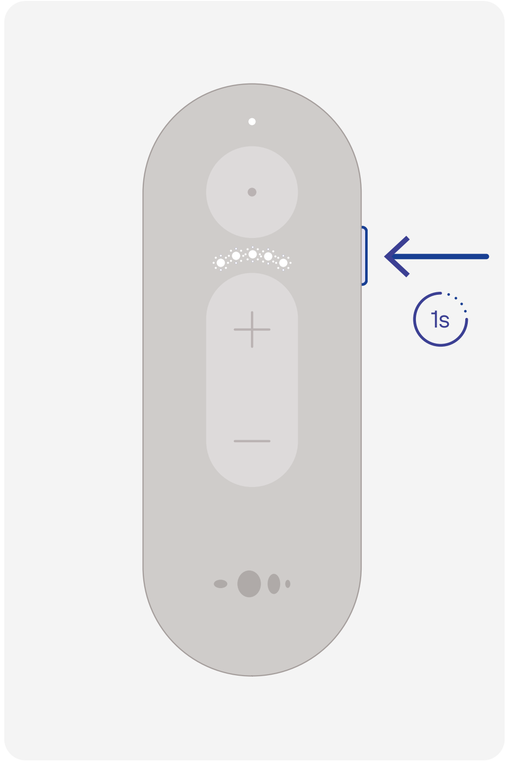


Cycle preset up
Press side-button and hold for 1 second.

Remote Plus user interface
Call

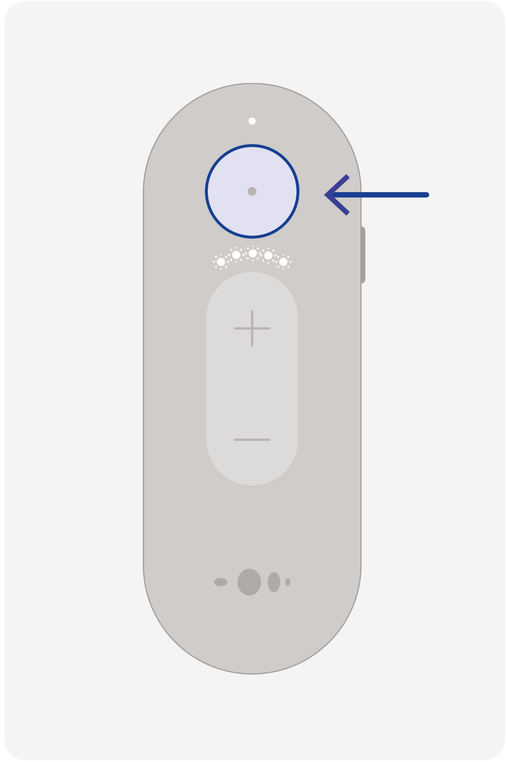

Answer call
Press top button

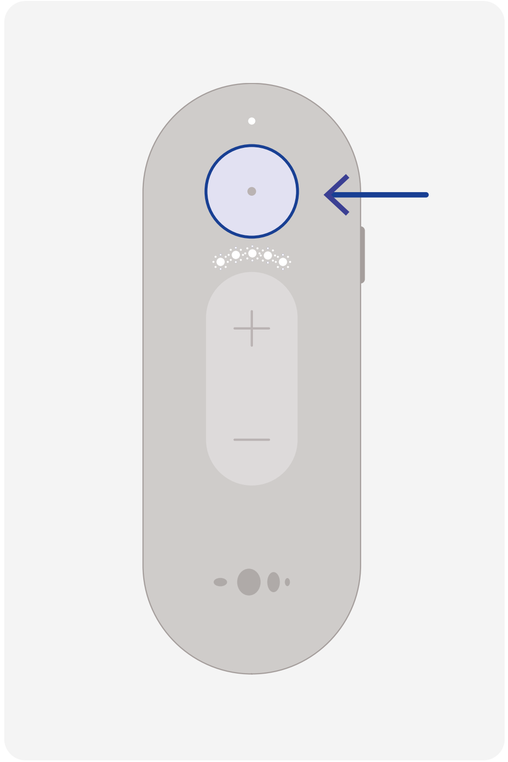

End call
Press top button

nuance audio

in tune with life

For the most updated Instructions for Use please follow this link

↗ View instruction for use